

Thai national guideline for nuclear medicine investigation in movement disorders: Nuclear medicine society of Thailand, the neurological society of Thailand, and Thai medical physicist society collaboration

Tawika Kaewchur^{1,2}, Benjapa khiewvan³, Wichana Chamroonrat⁴, Praween Lolekha⁵, Onanong Phokaewvarangkul⁶, Tanyaluck Thientunyakit³, Nantaporn Wongsurawat⁷, Peerapon Kiatkittikul⁸, Chanisa Chotipanich⁸, Wen-Sheng Huang⁹, Panya Pasawang¹⁰, Tanawat Sontrapornpol¹⁰, Nucharee Poon-iad³, Sasithorn Amnuaywattakorn⁴, Supatporn Tepmongkol^{11, 12, 13, 14*}

¹Division of Nuclear Medicine, Department of Radiology, Faculty of Medicine, Chiang Mai University Thailand

²PET/CT and Cyclotron Center, Center for Medicine Excellence, Faculty of Medicine, Chiang Mai University Thailand

³Division of Nuclear Medicine, Department of Radiology, Faculty of Medicine Siriraj Hospital, Mahidol University, Bangkok Thailand

⁴Division of Nuclear Medicine, Department of Diagnostic and Therapeutic Radiology, Faculty of Medicine Ramathibodi Hospital, Mahidol University, Bangkok, Thailand

⁵Division of Neurology, Department of Internal Medicine, Faculty of Medicine, Thammasat University, Pathumthani Thailand

⁶Department of Medicine, Faculty of Medicine, Chulalongkorn University, Bangkok, Thailand

⁷Division of Nuclear Medicine, Department of Radiology, Faculty of Medicine Siriraj Hospital, Mahidol University, Bangkok, Thailand

⁸National Cyclotron and PET Centre, Chulabhorn Hospital, Chulabhorn Royal Academy, Bangkok, Thailand

⁹Department of Nuclear Medicine, Show Chwan Memorial Hospital, Changhua, 500 Taiwan (R.O.C.)

¹⁰Division of Nuclear Medicine, Department of Radiology, King Chulalongkorn Memorial Hospital, Bangkok, Thailand

¹¹Division of Nuclear Medicine, Department of Radiology, Faculty of Medicine, Chulalongkorn University, Bangkok, Thailand

¹²Chulalongkorn University Biomedical Imaging Group (CUBIG), Faculty of Medicine, Chulalongkorn University, Bangkok, Thailand

¹³Chula Neuroscience Center, King Chulalongkorn Memorial Hospital, Thai Red Cross Society, Bangkok, Thailand

¹⁴Center of Excellence in Cognitive Impairment and Dementia, Faculty of Medicine, Chulalongkorn University, Bangkok, Thailand

ARTICLE INFO

Article type:
Others

Article history:

Received: 25 Oct 2023

Revised: 28 Nov 2023

Accepted: 6 Dec 2023

Keywords:

Guidance

Procedure Standard

Nuclear Neurology

Parkinson

ABSTRACT

Movement disorders are chronic neurological syndromes with both treatable and non-treatable causes. The top causes of movement disorders are Parkinson's disease and related disorders. Functional imaging investigations with Single Photon Emission Computed Tomography (SPECT) and Positron Emission Tomography (PET) images play vital roles in diagnosis and differential diagnosis to guide disease management. Since there have been new advanced imaging technologies and radiopharmaceuticals development, there is a need for up-to-date consensus guidelines. Thus, the Nuclear Medicine Society of Thailand, the Neurological Society of Thailand, and the Thai Medical Physicist Society collaborated to establish the guideline for Nuclear Medicine investigations in movement disorder for practical use in patient care. We have extensively reviewed the current practice guidelines from other related societies and good quality papers as well as our own experience in Nuclear Medicine practice in movement disorders. We also adjust for the most suitability for application in Thailand and other developing countries.

► Please cite this paper as:

Kaewchur T, khiewvan B, Chamroonrat W, Lolekha P, Phokaewvarangkul O, Thientunyakit T, Wongsurawat N, Kiatkittikul P, Chotipanich Ch, Huang W Sh, Pasawang P, Sontrapornpol T, Amnuaywattakorn S, Tepmongkol S. Thai national guideline for nuclear medicine investigation in movement disorders: Nuclear medicine society of Thailand, the neurological society of Thailand, and Thai medical physicist society collaboration. *Asia Ocean J Nucl Med Biol.* 2024; 12(2): doi:

* *Corresponding author:* Supatporn Tepmongkol. Division of Nuclear Medicine, Department of Radiology, Faculty of Medicine, Chulalongkorn University, 1873 Rama IV Road, Pathumwan, Bangkok, Thailand. 10330. Tel: +66 81 7743767; Email: supatporn.t@chula.ac.th

© 2024 *mums.ac.ir* All rights reserved.

This is an Open Access article distributed under the terms of the Creative Commons Attribution License (<http://creativecommons.org/licenses/by/3.0>), which permits unrestricted use, distribution, and reproduction in any medium, provided the original work is properly cited.

Introduction

Movement disorders are the clinical syndrome of neurological conditions affecting speed, fluency, quality, and ability to control movement. Movement disorders could either be an excess of movement or a paucity of voluntary and involuntary movements unrelated to weakness or spasticity. Movement disorders are often induced by pathological changes within the brain, especially in the basal ganglia and its connection circuits. The cause of movement disorders may be unknown, although hereditary and environmental factors may play a role in some cases. In this national guideline for nuclear medicine, we focus on the clinical use of functional imaging with nuclear medicine techniques in idiopathic Parkinson's disease (PD) and related disorders.

PD is one of the most recognizable movement disorders of the central nervous system, mainly affecting the motor system. The prevalence of PD is approximately 1% of the population above 60 years and increases with age (1, 2). The diagnosis of PD is based on a neurological examination with cardinal motor features: bradykinesia, rigidity, rest tremor, and postural instability (3). Since many other disorders can mimic signs and symptoms like those of PD called Parkinsonism, the differential diagnosis with other forms of Parkinsonism can be challenging, especially early in the disease, when signs and symptoms of other Parkinsonism are not fully shown. A prospective clinicopathological study showed a diagnostic accuracy of 65% within five years of disease onset (4). A clear response to dopamine replacement therapy in an individual with PD could help to reinforce that an accurate diagnosis has been established. However, acute challenge testing with levodopa should not be

used in the diagnosis of PD. Patients with suspected PD should be considered for a trial of chronic levodopa treatment (5). As the initial clinical diagnosis of PD is not fully precise, radiotracer neuroimaging could be an option for clinicians to improve diagnostic accuracy, which leads to better and more appropriate symptomatic management of patients, especially in the early stage of the disease (6). It has been well established that PD is associated with several non-motor symptoms. Some of these symptoms, such as hyposmia, constipation, depression, and REM sleep behavior disorder (RBD), might precede the onset of classical motor symptoms. Currently, there is no diagnostic test for PD.

Anatomical imaging, including Computed Tomography (CT) and Magnetic Resonance Imaging (MRI), is mainly used to exclude common secondary causes of Parkinsonism.

A recent revision of the diagnostic criteria of PD has gained attention on functional images (7).

Although they are not directly implemented in the clinical diagnostic criteria for PD, they are an optional assessment to evaluate atypical Parkinsonism. The scintigraphic evidence of cardiac sympathetic denervation on radiolabelled metaiodobenzylguanidine (MIBG) scan and normal presynaptic dopaminergic imaging on Positron Emission Tomography (PET) or Single Photon Emission Computerized Tomography (SPECT) are included in the supportive and absolute exclusion diagnostic criteria for PD, respectively (7). Furthermore, functional neuroimaging is considered a marker of the preclinical or prodromal PD stage (7-9). In addition, it could be used in the differential diagnosis of PD from other parkinsonian syndromes (Table 1) and monitoring the PD progression in the research study.

Table 1. Causes of Parkinsonism

Idiopathic Parkinson's disease (PD)	
Secondary parkinsonism	<ul style="list-style-type: none"> Drug-induced Vascular parkinsonism Normal-pressure hydrocephalus (NPH) Others; tumors, infection
Atypical parkinsonism	<ul style="list-style-type: none"> Progressive supranuclear palsy (PSP) Corticobasal degeneration (CBD) Multiple system atrophy (MSA) Dementia with Lewy bodies (DLB)
Heredodegenerative parkinsonism	

Presynaptic dopaminergic imaging has been used to demonstrate dopaminergic dysfunction underlying PD symptoms. Among these imaging modalities, ^{18}F -FDOPA PET and $^{99\text{m}}\text{Tc}$ -TRODAT-1 SPECT are currently available in Thailand. Fluorodopa (FDOPA) PET is used to measure aromatic amino acid decarboxylase (AADC) activity and reflect the dopaminergic nerve terminal integrity. At the same time, $^{99\text{m}}\text{Tc}$ -TRODAT-1 SPECT demonstrates the availability of presynaptic dopamine transporter (DAT).

However, the enzymatic / transporter targets visualized by FDOPA and TRODAT are often regulated by intrasynaptic dopamine concentration, and they can under- or over-estimate dopaminergic terminal densities. In PD, both techniques show a characteristic pattern of asymmetric reduction in striatal uptake that affects the posterior putamen contralateral to the clinically affected limb with a relatively preserved head of caudate (10-12). The reductions in tracer uptake also correlate with nigral dopamine cell counts and the degree of motor dysfunction, especially with bradykinesia. The findings mentioned above can also be found in atypical Parkinsonism. Thus, they cannot reliably differentiated among PD and others (13-15).

Regional cerebral glucose metabolism by fluorodeoxyglucose (^{18}F -FDG) PET reflects neuronal activity. Therefore, in the case of degenerative Parkinsonism, FDG PET might distinguish iPD patients from those with atypical Parkinsonism, such as multiple system atrophy (MSA), progressive supranuclear palsy (PSP), corticobasal degeneration (CBD), and dementia with Lewy bodies (DLB) (10, 11).

Also, the cortical glucose hypometabolism in parkinsonian patients might suggest cortical Lewy body dementia, such as Parkinson's disease dementia (PDD) and DLB.

Myocardial innervation scintigraphy with radio-labeled metaiodobenzylguanidine (MIBG) is a marker of the integrity and distribution of the postganglionic sympathetic nerve terminals.

Thus, in PD, the cardiac sympathetic denervation, as shown by the reduction or absence of MIBG uptake, can be used as a biomarker in the supportive criteria for diagnosing PD and differential diagnosis from other atypical Parkinsonism (12, 16).

To effectively use functional imaging with nuclear medicine, all members, including physicians, neurologists, radiologists, nuclear medicine professionals, and technicians, should have consensus information on these imaging techniques. In addition, functional imaging with nuclear medicine in movement disorders should be considered an aid to clinical diagnosis

in patients with uncertain diagnoses. However, there are limitations to accessing PET/CT services in some regions of Thailand. We thus propose alternative brain SPECT studies for each purpose.

Parkinsonism

Parkinsonism is a clinical syndrome characterized by the following cardinal features: rigidity, bradykinesia and rest tremor, typically asymmetrically affecting an upper extremity at onset (17). Postural instability is uncommon at the onset of PD but usually develops within several years of disease progression. Parkinsonism is a term used to describe movement symptoms associated with several conditions, including Idiopathic Parkinson's disease, atypical Parkinsonism (Parkinsonism Plus), and secondary Parkinsonism. Differences in abnormal protein are related to degenerative Parkinsonism and represented by various clinical parkinsonian syndromes such as synucleinopathies (PD, DLB, MSA) and tauopathies (PSP and CBD). The most reliable clinical diagnosis of PD is made with typical symptoms and signs (above) together with significant improvement with dopamine replacement medical therapy.

Idiopathic Parkinson's disease (PD)

In general, PD can be divided according to their manifestations into Prodromal and Clinical PD stages.

Prodromal Parkinson's disease

Prodromal-PD refers to the stage at which individuals have symptoms and signs that indicate a higher than average risk of developing PD (such as hyposmia, constipation, and REM sleep disorder), but do not manifest Parkinsonism sufficient to diagnose PD according to the diagnostic criteria (7).

Clinical Parkinson's disease

Clinical PD refers to the stage at which individuals present classical motor signs (bradykinesia, rigidity, and rest tremor) of PD and are diagnosed as PD according to the diagnosis criteria (3, 7, 8).

Parkinson's Disease Dementia (PDD)

Parkinson's disease dementia refers to the stage at which individuals with PD, who are diagnosed with PD according to established criteria, present with insidious onset and slow progression of cognitive impairment or dementia syndromes (18).

Atypical Parkinsonism

Atypical Parkinsonism previously known as a Parkinsonism-plus syndrome is a sporadic entity of neurodegenerative parkinsonian syndromes that could be distinguished from Parkinson's disease by a relative lack of response to dopamine replacement therapy and their (eventual) additional or "plus" features including, vertical supranuclear palsy (PSP), autonomic dysfunctions (MSA), dystonia or myoclonus or cerebral cortical signs (CBD), and early cognitive impairment (DLB) (19).

Therefore, it is necessary to identify red flag signs that should not be found in typical PD, and these signs support the diagnosis of atypical Parkinsonism.

Progressive Supranuclear Palsy (PSP)

In contrast to PD, PSP or Steele- Richardson-Olszewski syndrome presents a symmetrical akinetic-rigid syndrome. Rigidity in PSP usually predominantly involves the trunk and neck, causing early postural and gait instability with falls (20). Therefore, any elderly patient with truncal or neck rigidity with minimal limb involvement and impaired postural reflexes should be suspected of PSP. Vertical supranuclear gaze palsy is a pathognomonic sign of PSP, especially in a downward position (20).

Impairment of frontal executive functions is frequent in PSP, and dementia, often of a frontal type, is common in later diseases leading to memory impairment, expressive dysphasia, and limb apraxia.

Multiple System Atrophy (MSA)

MSA has two subtypes: predominant Parkinsonism (MSA-P, 60%) or predominant cerebellar symptoms (MSA-C, 40%). In addition to motor impairment, autonomic dysfunctions are a compulsory component of MSA-especially urinary incontinence, erectile dysfunction, or orthostatic hypotension (21).

MSA-P presents with a hypokinetic-rigid parkinsonian syndrome that tends to be more symmetrical and less responsive to levodopa than PD. In contrast, MSA-C usually presents with signs of cerebellar dysfunction including, limb and gait ataxia, scanning dysarthria, and nystagmus. The diagnostic criteria for MSA require at least one symptom of autonomic dysfunction (orthostatic hypotension, urinary incontinence, erectile dysfunction) after excluding alternative symptomatic causes (21).

Corticobasal Degeneration (CBD)

CBD is progressive atypical Parkinsonism that results from neurodegeneration of the basal ganglia and the cerebral cortex. The

combination of hypokinetic-rigid parkinsonian syndrome with additional signs, including dystonia, myoclonus, and cortical symptoms, is the characteristic feature of CBD, and is typically resistant to levodopa therapy (19).

Most patients eventually develop progressive immobility and cognitive impairment in a couple of years after diagnosis.

Dementia with Lewy Bodies (DLB)

Dementia with Lewy Bodies refers to a dementia syndrome with underlying cortical Lewy-related pathology that precedes parkinsonian motor signs or begins within one year from its onset (22). The definition of DLB is a progressive cognitive decline of sufficient magnitude to interfere with normal social or occupational functions or with usual daily activities. Prominent or persistent memory impairment may not necessarily occur early but is usually evident with progression (22). Deficits on attention, executive function tests, and visuo-perceptual ability may be especially prominent and occur early (22). In addition, four clinical core features are identified and usually occur and may persist throughout the course including, 1) fluctuating cognition with pronounced variations in attention and alertness, 2) recurrent visual hallucinations that are typically well formed and detailed, and 3) REM sleep behavior disorder, which may precede cognitive decline and 4) one or more spontaneous cardinal features of Parkinsonism.

Drug-Induced Parkinsonism (DIP)

Drug-induced Parkinsonism is the most common movement disorder induced by drugs that affect dopamine receptors (23). Since the clinical manifestations of DIP are very similar to those of Parkinson's disease (PD), patients with DIP are frequently misdiagnosed as having PD.

Typical antipsychotics, also known as neuroleptics, are the most common causes of DIP. In addition to antipsychotics, gastrointestinal (GI) motility drugs, calcium channel blockers (CCBs), and antiepileptic drugs have been found to induce DIP. Dopamine transporter (DAT) imaging may be used to diagnose various etiologies of Parkinsonism, including DIP (24).

Scans Without Evidence of Dopaminergic Deficit (SWEDDs)

SWEDDs are a heterogeneous group of patients who fulfill clinical diagnostic criteria for PD, but have normal dopaminergic functional imaging (25). Therefore, alternative diagnoses of SWEDDs have been suggested include essential tremor (ET), depression,

vascular or psychogenic Parkinsonism, dopa-responsive dystonia, and rare causes of supranigral Parkinsonism.

medical physicists and technologists as well as neurologists was established for guideline development. The involving panelists are listed in table 2. The guideline development process was shown in figure 1.

Guideline methodology and panelists

A panel of specialists in nuclear medicine,

Table 2. Panelist composition

Name	Field of expertise	Institution	Country
Supatporn Tepmongkol	Nuclear Medicine	Chulalongkorn University	Thailand
Tawika Kaewchur	Nuclear Medicine	Chiang Mai University	Thailand
Benjapa Khiewvan	Nuclear Medicine	Mahidol University	Thailand
Wichana Chamroonrat	Nuclear Medicine	Mahidol University	Thailand
Tanyalak Thientunyakit	Nuclear Medicine	Mahidol University	Thailand
Nantaporn Wongsurawat	Nuclear Medicine	Khon Kaen University	Thailand
Peerapon Kiatkittikul	Nuclear Medicine	National Cyclotron and PET Centre	Thailand
Chanisa Chotipanich	Nuclear Medicine	National Cyclotron and PET Centre	Thailand
Wen-Sheng Huang	Nuclear Medicine	Taipei Medical University Hospital	Taiwan
Panya Pasawan	Medical Physicist and Technologist	Chulalongkorn University	Thailand
Tanawat Sontrapornpol	Medical Physicist and Technologist	Chulalongkorn University	Thailand
Nucharee Poon-iad	Medical Physicist and Technologist	Mahidol University	Thailand
Sasithorn Amnuaywattakorn	Medical Physicist and Technologist	Mahidol University	Thailand
Praween Lolekha	Neurologist	Thammasat University	Thailand
Onanong Phokaewwarangkul	Neurologist	Chulalongkorn University	Thailand

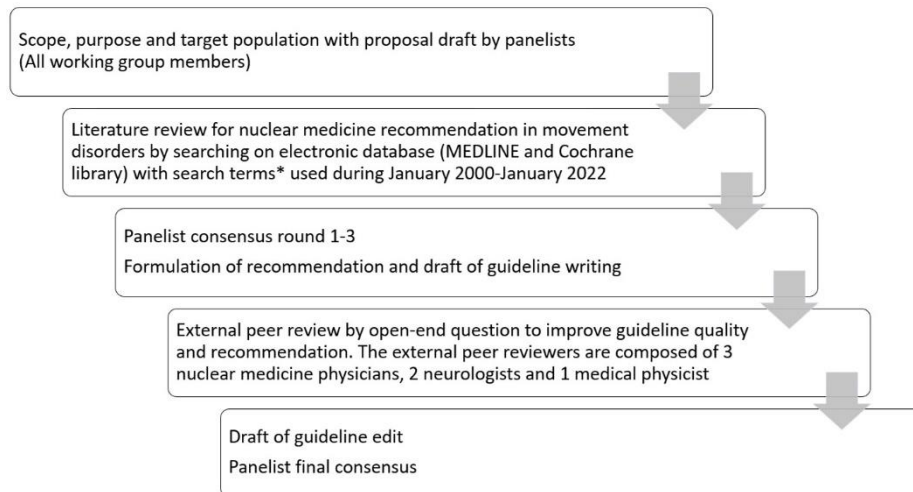


Figure 1. Guideline development process

Note :

1. Search terms are Parkinsonism, Parkinson’s disease, scintigraphy, SPECT, PET, diagnosis, differential diagnosis, and progression
2. All aforementioned search terms were used for keyword and mapping with medical subject headings in the hierarchy of central nervous system and movement disorders
3. Example of search formula ;
 Concept 1: Parkinsonism OR possible synonyms (eg. Parkinson plus syndrome, corticobasal degeneration)
 Concept 2: Scintigraphy OR SPECT OR PET
 Concept 3: Diagnosis OR differential diagnosis OR progression
 Then, searches these combinations to find all articles that include a search term from each concept: Concept 1 AND Concept 2 AND Concept 3

Objectives

1. To identify the role of nuclear medicine studies in patients with Parkinsonism.
2. To assist referrers in requesting the most appropriate procedure for Parkinsonism patients.

3. To identify scientific evidence that is useful to assist nuclear medicine professionals in recommending, performing, interpreting, and reporting the results of nuclear medicine investigations in Parkinsonism patients.

Indications and Contraindications

Indications for Nuclear Medicine Imaging in Parkinsonism with the Evidence Base and Flow of Investigation

The indication for nuclear medicine imaging in Parkinsonism was reviewed base on the established criteria for rating diagnostic accuracy studies and recommendations (American Academy of Neurology, 2017). The evidence-based reviews have used a standard method using literature searches performed using electronic databases (Medline, Cochrane Library) with search keywords: parkinsonism, Parkinson's disease, scintigraphy, SPECT, PET, diagnosis, progression, differential diagnosis, from January 2006 to June 2022. The included

studies have to be human-based with full English articles. Each study was rated by 2 independent panelists (P.L. and O.P.). The overall recommendation for clinical practice was then assessed and classified as A: established effective, B: probably effective, C: possible effective, and U: data inadequate.

The evidence of nuclear medicine investigation in Parkinsonism is shown in table 3 and the flow of imaging investigation in movement disorder is shown in Figure 2.

Contraindications

- Pregnancy
- Breastfeeding

Table 3. Level of Recommendation for Nuclear Medicine Imaging in Parkinsonism

	Level of Recommendation
Support a diagnosis of PD	
FDOPA PET	B
DAT SPECT/PET	A
FDG PET	U
Cardiac MIBG	U
Assessment of the progression of PD	
FDOPA PET	B
DAT SPECT/PET	B
FDG PET	U
Cardiac MIBG	U
Differentiate PD from atypical parkinsonism	
FDOPA PET	U
DAT SPECT/PET	U
FDG PET	C
Cardiac MIBG	U

A: established effective, B: probably effective, C: possibly effective, and U: data inadequate

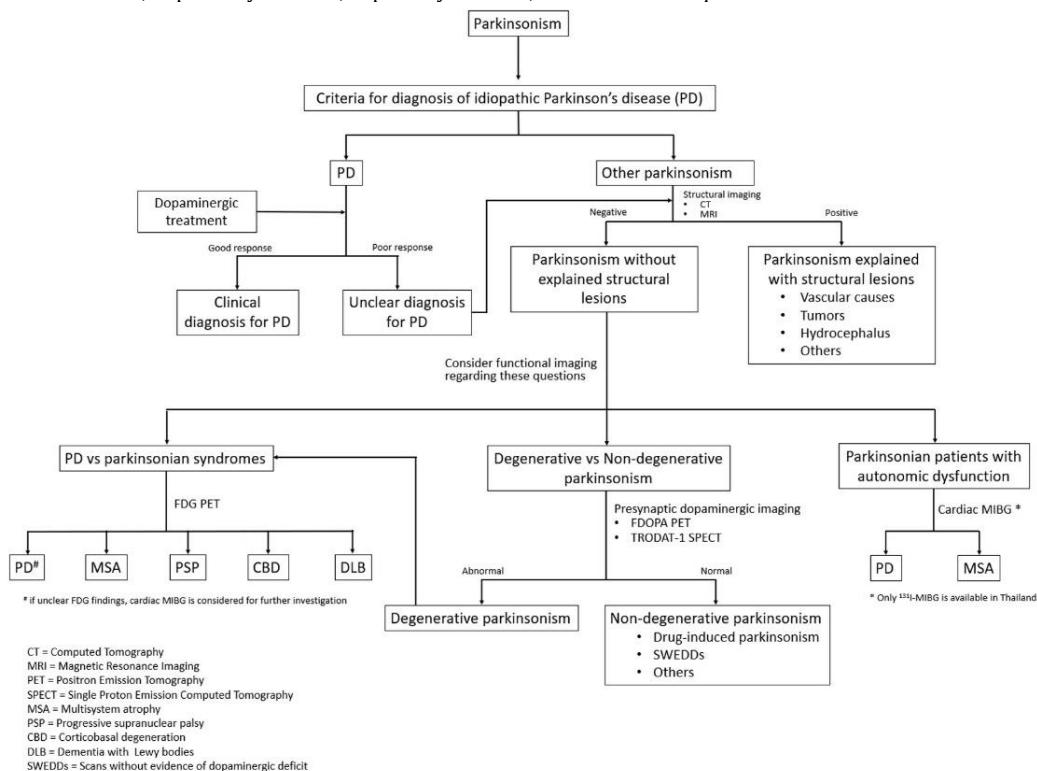


Figure 2. Flow of Investigation in Movement Disorder

Procedure

Patient Preparation

I. TRODAT-1 SPECT

Patient Preparation and Precautions (26)

Check for medications on table 4 that may alter tracer binding and need to withdraw (if possible) for at least five half-lives.

Table 4. Lists of Medication that May Alter TRODAT-1 Uptake

Medication	Effect of the medication
Cocaine	Severely decrease binding to DaT
Amphetamines	
<ul style="list-style-type: none"> d-amphetamine methamphetamine methylphenidate 	Decrease striatal binding
Central nervous system stimulants	
<ul style="list-style-type: none"> ephedrine phentermine 	Decrease striatal binding
*particularly when used as tablets	
Modafinil	Decrease striatal binding
Some antidepressants (mazindol, bupropion, radafaxine)	Decrease striatal binding
Adrenergic agonists	
<ul style="list-style-type: none"> phenylephrine norepinephrine 	Increase striatal binding
Anticholinergic drugs	
<ul style="list-style-type: none"> benztropine other anticholinergics 	Decrease striatal binding ratios Might increase these ratios, which will likely not affect visual assessments
Opioids (only fentanyl)	Decrease striatal binding
Bupropion, fentanyl, and some anesthetics (ketamine, phencyclidine, and isoflurane)	May decrease binding to DaT
Selective serotonin reuptake inhibitors	May increase binding to DaT somewhat but should not interfere with visual interpretation
Cholinesterase inhibitors and neuroleptics	Probably do not interfere significantly with binding to DaT
Anti-Parkinsonian drugs (e.g., L-dihydroxyphenylalanine, dopamine agonists, monoamine oxidase B inhibitors, N-methyl-D-aspartate receptor blockers, and catechol-O-methyltransferase inhibitors in standard dosages)	Do not interfere with binding to DaT to any significant degree

II. Cardiac MIBG Imaging

Patient Preparation and Precautions (27)

1. Withdrawal of medications listed in table 5, which may alter tracer uptake, is advised. In the case of hypertension, the recommended anti-hypertensive medications are shown in Table 6.
2. Thyroid Blockade
The thyroid blockade agent must be given one day before tracer injection and continue for one to two days for ¹²³I-MIBG or two to three days for ¹³¹I-MIBG (27). The recommendations for

thyroid blockade are shown in table 7.

On the Acquisition Day

1. Evaluate baseline clinical, vital signs, and general physical examination.
2. Slowly inject MIBG and observe pallor, vomiting, or fainting.
3. Encourage patients to drink a large volume of fluid and void before the study .
Monitor vital signs, and observe the adverse reactions of MIBG (tachycardia, pallor, vomiting, abdominal pain).

Table 5. Medication and drugs of abuse interfering I-131 MIBG uptake

Medication	Approved name	Effect of the medication	Withdraw period recommendation in tumor
I. Cardiovascular and sympathomimetic drugs			
Antiarrhythmics for ventricular arrhythmias	Amiodarone	Inhibition uptake and depletion	Not practical to withdraw
Combined α/β -blocker	Labetalol*	Inhibition uptake and depletion	72 hours 21 days*
Adrenergic neuron blocker	Bretylium*, Guanethidine*, Reserpine*	Depletion and transport inhibition	48 hours 14 days*
α -Blocker	Phenoxybenzamine (intravenous doses only)	Unknown mechanisms	15 days
Calcium channel blockers	Diltiazem, Nifedipine*, Nimodipine	Increased uptake and retention	24 hours 14 days*
	Amlodipine*, Felodipine, Isradipine, Lacidipine, Lercanidipine, Nicardipine*, Nisoldipine, Verapamil	Increased uptake and retention	48 hours 14 days*
Inotropic sympathomimetics	Dobutamine, Dopamine, Dopexamine	Depletion of granules	24 hours
Vasoconstrictor sympathomimetics	Ephedrine	Uptake inhibition	24 hours
	Metaraminol, Norepinephrine, Phenylephrine	Depletion of granules	24 hours
β 2 stimulants (sympathomimetics)	Salbutamol, Terbutaline, Eformoterol, Bambuterol, Fenoterol and Salmeterol	Depletion of granules	24 hours
Other adrenoceptor stimulants	Orciprenaline	Depletion of granules	24 hours
Systemic and local nasal decongestants, compound cough and cold preparations	Ephedrine	Uptake inhibition	24 hours
	Xylometazoline, Oxymetazoline	Depletion of granules	24 hours
	Pseudoephedrine, Phenylephrine	Depletion of granules	48 hours
Sympathomimetics	Brimonidine, Dipivefrine, Phenylpropanolamine*, Ephedrine*, Pseudoephedrine*, Phenylephrine*, Amphetamine*, Dopamine*, Isoproterenol*, Salbutamol*, Terbutaline*, Phenoterol*, Xylometazoline*	Depletion of granules	48 hours 7-14 days*
II. Neurological drugs			
Antipsychotics (neuroleptics)	Chlorpromazine*, Perphenazine, Prochlorperazine, Promazine, Thioridazine	Uptake inhibition	24 hours 21-28 days*
	Fluphenazine*	Uptake inhibition	24 hours, or 1 month for depot 21-28 days*
	Benperidol, Pericyazine, Sulpiride, Trifluoperazine, Quetiapine	Uptake inhibition	48 hours
	Droperidol*, Flupentixol, Haloperidol*, Zuclopenthixol	Uptake inhibition	48 hours, or 1 month for depot 21-28 days*
	Levomepromazine, Pimozide, Amisulpride	Uptake inhibition	72 hours
	Zotepine	Uptake inhibition	5 days
	Risperidone	Uptake inhibition	5 days or 1 month for depot
	Clozapine	Uptake inhibition	7 days
	Olanzapine	Uptake inhibition	7-10 days
	Sertindole	Uptake inhibition	15 days
	Pipotiazine	Uptake inhibition	1 month for depot
	Lozapine*	Uptake inhibition	21-28 days*
	Sedating antihistamines	Promethazine*	Uptake inhibition
Opioid analgesics*	Tramadol*	Uptake inhibition	24 hours 7-14 days*
Tricyclic antidepressants	Clomipramine, Dosulepin (dothiepin), Doxepin*, Imipramine*, Nortriptyline	Uptake inhibition	24 hours 7-14 days*
	Amitriptyline*, Amoxapine*, Lofepramine, Trimipramine	Uptake inhibition	48 hours 7-14 days*
Tricyclic-related antidepressants	Maprotiline*, Mianserin, Trazolone*, Venlafaxine	Uptake inhibition	48 hours 21-28 days*
	Reboxetine	Uptake inhibition	72 hours
	Mirtazepine	Uptake inhibition	8 days
CNS stimulants	Cocaine*, Caffeine	Uptake inhibition	24 hours 7-14 days*
	Amphetamines (e.g. dexamfetamine), Methylphenidate	Depletion of granules	48 hours
	Modafinil	Unknown mechanisms	72 hours
	Atomoxetine	Uptake inhibition	5 days

NOTE:

Modified from;

EANM guideline >>> 131I-/123I-Metaiodobenzylguanidine (MIBG) scintigraphy: procedure guidelines for tumor imaging; Published online: 20 July 2010, EANM 20103

* Cardiac MIBG guideline and study review⁴. Withdrawal period is preferred using the time interval for cardiac MIBG study. However, if it cannot be applied, using tumor recommendation withdrawal period is also acceptable

Table 6. Recommended Antihypertensive Drugs for 131I-MIBG Study

Medication	Example of Medications
α-blocker	Selective α1-blockers: Prazosin, Terazosin, Doxazosin Phenoxybenzamine
β-blocker	Propranolol

Table 7. Thyroid Blockade in Adults

Potassium iodate	170 mg
Potassium iodide	130 mg
Lugol's 1%	1 drop/kg to a maximum of 40 drops (20 drops twice a day)
Potassium perchlorate	400 mg

III. FDG PET

Patient Preparation and Precautions (28-30)

1. Patients should fast for at least four hours to allow optimal cerebral FDG uptake.
2. Advise the patients to avoid excessive stimulants which may affect brain metabolism at least 24 hours.
 - Food or drinks such as caffeine, cola, alcohol, and energy drinks.
 - Smoking and excessive exercise.
3. Control blood glucose in diabetic patients
4. Recommended medications to be withdrawn are shown in table 8.

On Acquisition Day

1. Blood glucose must be checked in every case. The optimal blood glucose level should not exceed 160 mg/dl (30). However, high blood glucose level exceeding 160 mg/dl is not a contraindication but needs a caution to interpret due to elevated brain background activity.

2. Adequate-hydration is advised.
3. Check the patient's ability to lie still during the whole acquisition time.
4. Intravenous catheterization should be placed at least ten minutes before injection.
5. Place the patient in a quiet, dimly lit room with open eyes and unplugged ears with minimal to no movement and no external stimulus for at least 30 minutes before and after tracer injection.
6. The optimal uptake time is 30-60 minutes.
7. Instruct the patient to void before scanning for patient comfort and to minimize radiation exposure to the bladder.
8. In uncooperative patients, sedation (28) can be applied and should be given just before image acquisition.
9. Continually monitor vital signs, including pulse oximetry, during the entire scanning procedure, especially in those receiving sedation. Emergency backup should be available.

Table 8. Medications interfering FDG uptake

Medication	Duration of withdrawal
Amphetamine	4-6 half-life
Cocaine	4-6 half-life
Anesthetics	
Propofol	4-6 half-life
Isoflurane	
Barbiturate	
Benzodiazepine	4-6 half-life
Corticosteroid	4-6 half-life
Neuroleptics (anti-psychotics)	
Haloperidol	4-6 half-life
Chlorpromazine	
Cholinesterase inhibitor	
Donepezil	4-6 half-life
Rivastigmine	

IV. FDOPA PET

Patient Preparation and Precautions (31)

1. Patients should fast for at least two to six hours. Do not take any amino-acid-containing foods six hours before the procedure.
2. Medication withdrawal: stop COMT inhibitors for five half-lives before scanning. Other drugs used for PD can be continued.

On Acquisition Day

1. Pre-medication: Pure carbidopa 2 mg/kg (max 150 mg) orally, 90 minutes prior to ¹⁸F-FDOPA injection is optional (32).

2. Pre-hydration: drink up to one liter of water before intravenous administering ¹⁸F-FDOPA (200 MBq) to decrease the ¹⁸F-FDOPA urinary concentration.
3. Conscious sedation (28) (e.g., a short-acting benzodiazepine such as intravenous midazolam) can be applied in uncooperative or uncontrolled-symptom patients. Sedation should be given at least 20 minutes after radiopharmaceutical administration.
4. Continually monitor vital signs, including pulse oximetry, during the entire scanning

procedure, especially in those receiving sedation. Emergency backup should be available.

Radiopharmaceuticals and radiation dosimetry

Radiopharmaceutical dosage is shown in Table 9 (28, 33-36), and their characteristics are shown on table 10.

Radiation Dosimetry

Radiation dosimetry in SPECT and PET studies is shown in Table 11 (29, 33, 35, 37).

SPECT Imaging Acquisition and Data Processing

TRODAT-1 Brain SPECT

The patient positioning, image acquisition and data processing of TRODAT-1 brain SPECT are shown in Table 12 (38).

Table 9. Radiopharmaceutical Dosages

Radiopharmaceuticals	Dosage
^{99m} Tc-TRODAT-1	740 MBq (20 mCi)
¹³¹ I-MIBG	111 MBq (3 mCi)
¹⁸ F-FDOPA	3.5 MBq/kg (0.09 mCi/kg)
¹⁸ F-FDG	185-370 MBq (5-10 mCi)

Note: Radiopharmaceuticals dose in this table is general recommendation. Dose adjustment should be considered for local regulatory approval as well as types of scanners used

Table 10. Radiopharmaceuticals Characteristics

Agent	Peak brain activity(min)	Brain uptake (% of injected dose)	Brain/Bg Ratio	GM: WM ratio	Imaging time	Brain washout
TRODAT-1	30-40 (striatum in baboon*)	<2	2:1 (basal ganglia: Occipital) 1.5-2.2:1 (Striatum: Cerebellum)	N/A	3.5-5 hours	Slow washout from the striatum after 90 mins
FDOPA	30-50	6	Higher	2.5-4.1:1	>30 min	none

Agent	Peak cardiac activity (min)	Cardiac uptake (%)	Early heart /mediastinum	Delayed heart /mediastinum ratio (4 hrs)	Imaging time	Cardiac washout
MIBG	120	< 1% at less than 24 hours	2.41 + 0.26	2.66 + 0.47	15-30 minutes for early 3-4 hours for delayed images	Yes

Table 11. Radiation Dosimetry in SPECT and PET Studies

Radiopharmaceuticals	Dose in MBq unit (mCi)	Organ receiving the highest dose mGy/MBq	Effective dose mSv/MBq
^{99m} Tc-TRODAT-1	740 (20)	0.047 (liver)	0.015
¹³¹ I-MIBG	44.4 – 81.4 (1.2 – 2.2)	0.83 (liver)	0.14
¹⁸ F-DOPA	200 (5.4)	0.15 (Bladder wall)	0.0199
¹⁸ F-FDG	185-740 (5-20)	0.13 (bladder wall)	0.019

Note: Low-dose CT scan is used for attenuation correction. The effective dose of a low dose CT scan is about 20 μSv and 220-450 μSv in a high dose or diagnostic CT

Table 12. TRODAT-1 SPECT Acquisition and Data Processing

Patient Positioning	
Supine and arms down position on the scanner table	
Immobilization to maintain head position and reduce movement	
Patient's head in the center of the field of view	
The head held straight (not tilted left or right) with the line from ear to the eye (the canthomeatal line) perpendicular to the detector surface, if possible	
Use the smallest radius of rotation as possible or automated contour setting from patient to detector for imaging to ensure maximum image resolution	
Data Acquisition Parameters	
Uptake time	4 hours
Instrument	SPECT
Collimator	Fan beam or Parallel hole (LEHR/LEUHR)*
Energy setting	140 keV, 15-20% energy window
Nuclide	^{99m} Tc
Matrix size	≥ 128×128
Zoom	To gain at least pixel size smaller than 1/3 to 1/2 of the expected resolution
Scan mode	Step and shoot or continuous
Rotation per view	≤ 3° (total of 360° rotation)
Time per view	About 15-30 sec/projection (minimum total count 5 x 10 ⁶ counts)
Scatter correction	Optional, but recommended
Reconstruction	3D-OSEM/FBP**
Slice thickness	Possible 3-5 mm (for maximal pixel resolution)
If SPECT/CT (CT for attenuation correction)	
CT voltage	120-140 kV
CT current	30-80 mA (depends on vender's recommendation)
Slice thickness	Same as SPECT slice thickness for AC
Data Processing	
Review projection data in cine mode and sonogram for an initial determination of image quality, patient motion, and artifacts	
Reconstruct by OSEM or FBP algorithm	
Select resolution recovery filter, i.e., Metz filter or type of low pass filter, i.e., Butterworth, Hamming, or Hanning	
Optimize reconstructing parameters, i.e., cut-off, order, iterations, and subset, depending on the injected activity, camera sensitivity, and resolution	
Apply attenuation correction using either calculated (e.g., Chang's method) or measured (e.g., Gadolinium source or CT scan) attenuation	
Apply the scatter correction (optional) to improve the image signal-to-noise ratio with various methods. The most popular one is the triple-energy window correction.	

Note:

*Low Energy High Resolution (LEHR) or Low Energy Ultra High Resolution (LEUHR)

** Ordered Subsets Expectation Maximization (OSEM) or Filtered Back Projection (FBP)

Cardiac MIBG Planar Imaging

The patient positioning, image acquisition and data processing of cardiac MIBG planar image are shown in Table 13.

Brain PET Imaging Acquisition and Data Processing

The patient positioning, image acquisition and data processing of FDG brain and FDOPA brain PET are shown in Table 14.

Table 13. Cardiac MIBG Planar Image Acquisition and Data Processing

Patient Positioning	
Comfortably supine and arms down position on the imaging table	
Immobilization to maintain chest position and reduce movement	
The patient's heart is in the center of the field of view, covering the mediastinum and whole lung	
Data Acquisition Parameters	
Uptake time	Early image at 15-30 minutes post tracer injection Delayed image at 3- 4 hours post tracer injection
Instrument	SPECT
Collimator	High Energy
Energy setting	364 keV, 15-20% energy window
Nuclide	¹³¹ I
Matrix	256×256
Mode	Static
View	Anterior
Zoom	1.67 - 2 (< FOV)
Scan time	5 - 10 min
Data Processing	
Display, create and record average count region of interest (ROI) cover myocardium and upper mediastinum in anterior view of early and delayed images. Then, calculate the heart to mediastinum ratio (H/M ratio, as the quotient [average count in the myocardium]/ [average count in the upper mediastinum]) at early and delayed images.	

Table 14. FDG Brain and FDOPA Brain PET Acquisition and Data Processing

Patient Positioning	
Comfortably supine and arms down position on the imaging table	
Place the patient's head in a holder using foam or folded sheet so that the canthomeatal line is in a vertical orientation	
Immobilization to maintain head position and reduce movement	
Patient's head in the center of the field of view	
Data Acquisition Parameters	
Uptake time	30-60 minutes for FDG PET 90 minutes for FDOPA PET
CT	
Scan type	Helical
Rotation time	0.75-1 section
Matrix	512
Slice thickness	3-5 mm
Slice increments	Continuous
Pitch	≤ 1
CT voltage	120-140 kV
CT current time	≤ 50 mAs (Low dose CT) ≤ 250 mAs (Diagnostic CT)
PET	
Energy setting	511 keV, 15-20% energy window
Nuclide	¹⁸ F
Mode	2D, 3D in list mode
Scan direction	Towards head
Scan duration	5-15 minutes for FDG PET (vary depending on vendor's recommendation) 30 minutes for static FDOPA PET image (may vary depending on vendor's recommendation) Dynamic FDOPA PET image is optional with 5 minutes/ frame, a total of 18 frames
CT Attenuation correction	Yes
Scatter correction	Yes
Reconstruction	Iterative (OSEM) or equivalent
Data Processing	
Preview image for patient motion and ensure PET and CT images are matched before performing CT attenuation correction	
Images are reconstructed in the transaxial plane of at least 128X128 matrix size, zoom 2	
The typical pixel size is 2-4 mm	
Depending on the resolution of the PET system, a final image resolution may vary between 2.5-10 mm. FWHM. This FWHM typically yields adequate image resolution and signal-to-noise ratios	
The reconstruction parameters can be varied. Please refer to the manufacturer's recommendations for the best choices of iterations, subsets, and smoothness.	

Interpretation

Cardiac MIBG

General Considerations

The radiolabeling of MIBG with iodine-123 or iodine-131 enables the visualization and quantification of cardiac postganglionic sympathetic innervation in vivo. Cardiac MIBG scintigraphy has demonstrated its capability to differentiate between PD and atypical Parkinsonism, essential tremor (ET), and normal control (33, 39-43). Furthermore, it can also differentiate between dementia with Lewy

bodies and other types of dementia (43). In Thailand, only the ¹³¹I isotope is available. Thus, ¹³¹I-MIBG is used in this guideline.

Image Data Display

The anterior view of chest planar imaging is displayed in grayscale at 15-30 minutes (early scan) and 3-4 hours (delayed scan). The ratio of ¹³¹I-MIBG uptake in the heart to the mediastinum (H/M ratio) is quantified by drawing the region of interest (ROI) around the heart and the upper mediastinum, as shown in Figure 3.

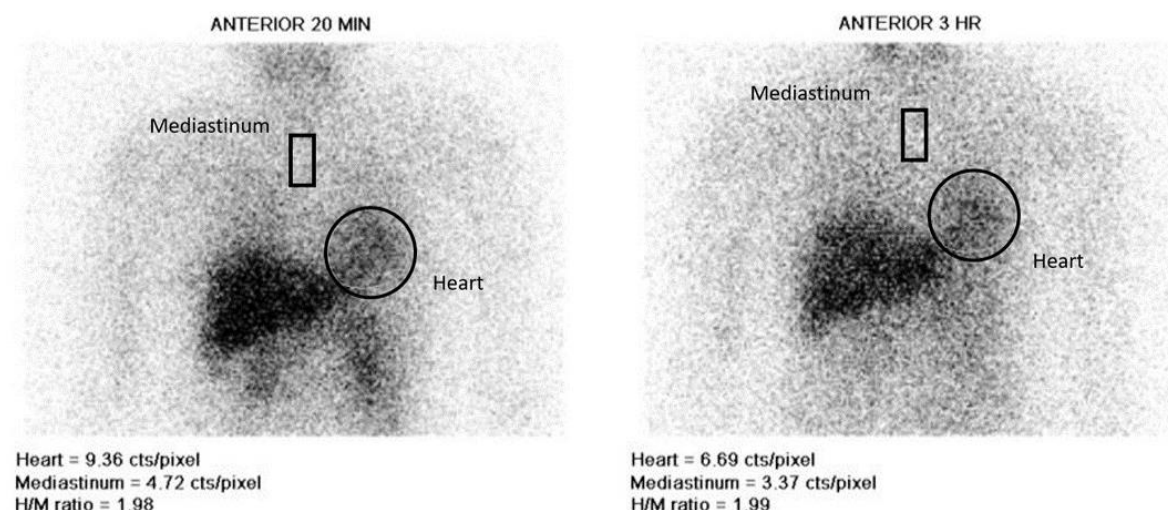


Figure 3. I-131 MIBG anterior view of the chest at 20 minutes and three hours. ROIs were drawn around the heart and mediastinum to calculate the H/M ratio

Image Interpretation

Semi-quantitative analysis using heart/mediastinum (H/M) ratio is often used for cardiac ¹³¹I-MIBG analysis. The mediastinum is often used as a reference because it contains few sympathetic nerves. Therefore, average counts per pixel in the heart and mediastinum are used for calculation. The washout ratio (WOR) of MIBG is another parameter used. It indicates the activity tone of the presynaptic sympathetic nerves, which is significantly higher in PD. The WOR is calculated as follows (33).

$$WOR = \frac{[Heart(early) - Mediastinum(early)] - [Heart(late) - Mediastinum(late)]}{[Heart(early) - Mediastinum(early)]}$$

Note: early = 15 minutes and late = 4 hours after ¹³¹I-MIBG administration

The H/M ratios at early (15 minutes) and late (4 hours) time points and WOR of ¹³¹I-MIBG are shown in Table 15 (33).

H/M ratio can also differentiate dementia with Lewy bodies (DLB) and Alzheimer’s disease (AD).

Since the Lewy bodies also causes DLB, myocardial uptake is also diminished. A study using ¹³¹I-MIBG (43) demonstrated that the H/M ratio at early (30 minutes) and late (4 hours) time points were lower in DLB than in AD, as shown in Table 16.

Table 15. Early (15 minutes) and late (4 hours) H/M ratio with standard deviation (SD) of 131I-MIBG in HC, ET, MSA, and PD

Condition	Early H/M ratio	Late H/M ratio	WOR
HC	2.41 ± 0.26	2.66 ± 0.47	-7.87 ± 32.65
ET	2.34 ± 0.34	2.46 ± 0.51	-4.70 ± 30.75
MSA	1.97 ± 0.36	2.08 ± 0.57	9.96 ± 30.88
PD	1.65 ± 0.36	1.50 ± 0.43	32.58 ± 22.99

Note: HC = healthy control

Table 16. Early (30 minutes) and late (4 hours) H/M ratio range of 131I-MIBG in Dementia with Lewy body (DLB) and Alzheimer’s disease (AD)

Condition	Early H/M ratio	Late H/M ratio
DLB	1.37-1.49	1.12-1.54
AD	1.72-1.79	1.78-1.80

TRODAT-1 brain SPECT

General Considerations

TRODAT SPECT represents the dopamine transporter. The normal study shows symmetrical radiotracer uptake in both caudate nuclei and putamina with very low or absent background activity. The head of caudate and putamen show the highest activity, as seen in the “comma” shape for all ages. However, the degree of striatal binding may decrease with normal aging by about 5-7% per decade (44).

However, the pattern for visual interpretation remains unchanged. TRODAT SPECT imaging is not recommended to distinguish between PD and atypical Parkinsonism and among atypical Parkinsonism.

Image Data Display

TRODAT-1 brain SPECT was displayed in the axial view in gray and color scale using the level of anterior commissure-posterior commissure (AC-PC) line. Raw data images must be checked

for movement correction. The regions of interest (ROIs) were bilateral striata, bilateral caudate nuclei, and bilateral putamina.

Image Interpretation

1. Visual Inspection (table 17 and figure 4)
2. Semi-Quantitative Analysis

Semi-quantitative analysis helps to determine disease accurately by calculating the specific uptake ratio (SUR) and asymmetrical percentage of bilateral striata (Figure 5). The percentage of asymmetry is significant if more than 10%, indicating early Parkinson’s disease (45). There is no consistently accepted cut-off value for SUR since the value can be influenced

by various factors including reconstructed image resolution, attenuation correction, and the size and the placement of regions of interest. The three or more transverse consecutive slices with the highest uptake at striatum were summed up and analyzed. The regions of interest were bilateral striata, bilateral caudate nuclei, and bilateral putamen. The size and shape of ROI should be fixed for each individual.

$$\text{Specific uptake ratio (SUR)} = \frac{(\text{counts at region of interest} - \text{occipital counts})}{\text{occipital counts}}$$

(46, 47)

Region of interest = caudate head, putamen, striatum

$$\text{Asymmetrical ratio of bilateral striata} = \frac{(\text{SUR right striatum} - \text{SUR left striatum}) \times 100\%}{(\text{SUR right striatum} + \text{SUR left striatum})} \quad 41$$

Table 17. TRODAT-1 Visual Interpretation

Visual score and findings	Possible disease
<p>Normal</p> <p>Class 0 = normal bilateral uptake of striatum</p>	<p>Normal</p> <p>ET</p> <p>drug-induced parkinsonism</p> <p>vascular parkinsonism</p> <p>psychogenic parkinsonism</p> <p>SWEDDs</p> <p>AD</p>
<p>Abnormal</p> <p>Class 1 = normal caudate uptake, loss putamen uptake > 50% on one side</p> <p>Class 2 = normal caudate uptake, loss putamen uptake > 50% both sides</p> <p>OR decreased caudate and putamen uptake about 50% on both sides</p> <p>Class 3 = loss caudate uptake > 50% and no putamen uptake</p> <p>Class 4 = finding between class 3 and 5</p> <p>Class 5 = no striatal uptake both sides</p>	<p>Early PD</p> <p>Early PD</p> <p>Early PD</p> <p>Advanced PD</p> <p>Advanced PD</p> <p>Advanced PD</p>

Note: Reduction of TRODAT-1 uptake in the striatum can be found in PD, DLB, atypical Parkinsonism such as PSP, MSA, and CBD

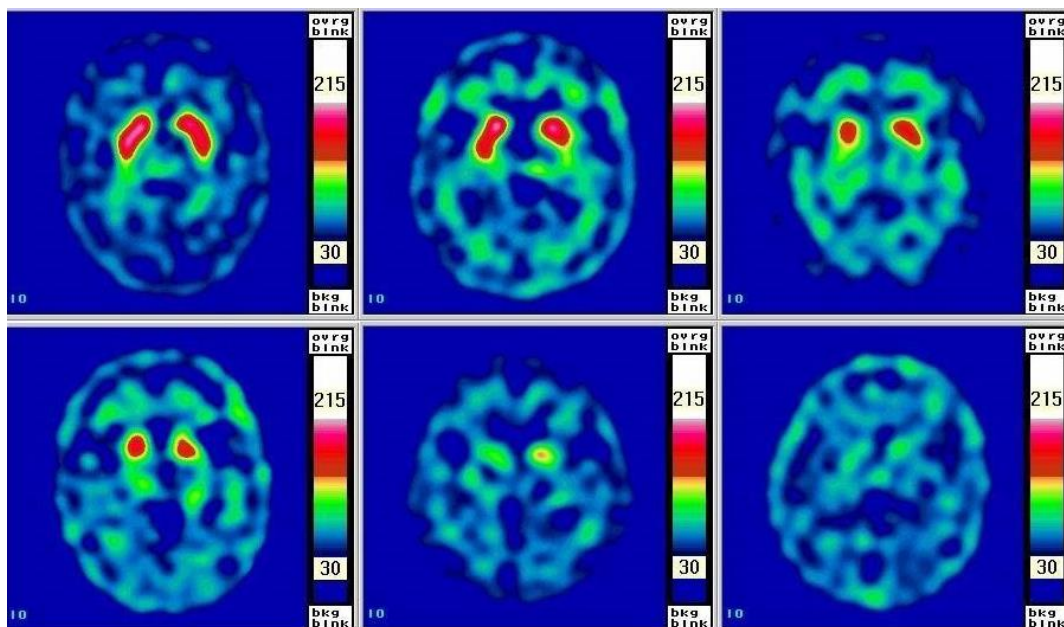


Figure 4. Normal to abnormal TRODAT-1 patterns in H-Y STAGE 0-5
 Reprinted by permission from: [Springer] [Eur J Nucl Med Mol Imaging] [Usefulness of brain ^{99m}Tc-TRODAT-1 SPET for the evaluation of Parkinson’s disease. Huang.WS, Lee.MS, Lin. JC et al.] [COPYRIGHT] (2004)

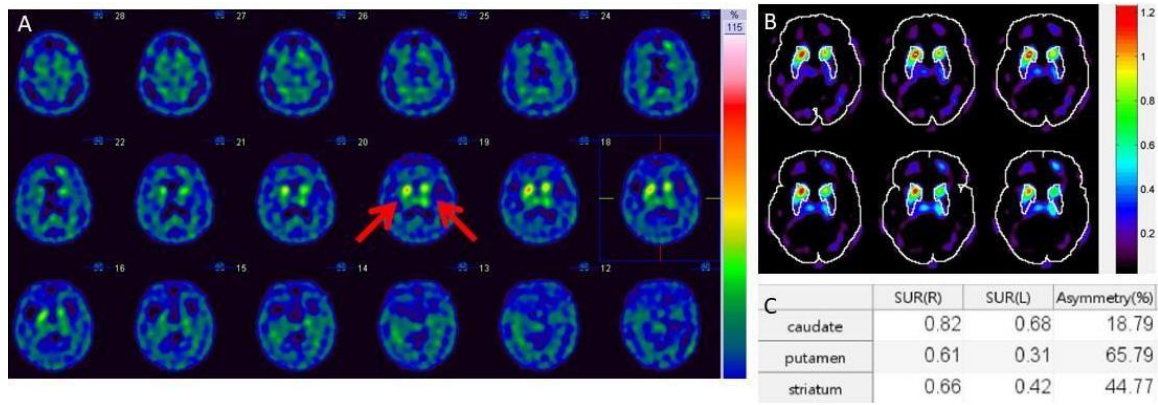


Figure 5. The TRODAT-1 SPECT in Parkinson’s disease with H-Y stage II. **A:** the TRODAT-1 SPECT showed moderately asymmetrical decreased DAT activity in bilateral putamina, and caudate nuclei. **B** and **C** showed a 65.79% asymmetry of putamen uptake. SUR of left putamen is 0.31, while SUR of right putamen is 0.61

FDG Brain PET

General Considerations

Brain glucose metabolic rate is considered to be a general marker of neuronal function. Thus, PET/CT imaging by ¹⁸F-FDG reflects regional neuronal function and provides the disease-specific patterns of resting regional glucose metabolic change accompanying movement disorders.

Image Data Display

1. Axial image in AC-PC (Anterior commissural-Posterior commissural plane) orientation as well as coronal and sagittal views.
2. Image co-registration with MRI image is an option for accurate anatomical localization.

Image Interpretation

Summarized findings of each movement disorder abnormality are shown in Table 18 and Figure 6(2, 48, 53).

Table 18. FDG Brain PET Findings in Parkinson’s disease and Parkinsonism

Disease	Findings
Parkinson disease	Increased FDG uptake: thalamus, striatum (or normal), sensorimotor cortex, cerebellum Decreased FDG uptake: temporal, parietal, possible occipital, and frontal cortices in demented PD (associated with future cognitive impairment)
MSA	Decreased FDG uptake MSA-P: striatum MSA-C: cerebellum
PSP	Decreased FDG uptake: dorsolateral frontal cortex, medial frontal cortex, thalamus, caudate nucleus, and upper brain stem
CBD	Decreased FDG uptake: asymmetrical at the striatum, medial, and lateral frontoparietal cortices (contralateral to the most affected body side)
DLB	Decreased FDG uptake: parietotemporal association cortices, occipital association cortex, and primary visual cortex (medial occipital lobe), bilaterally. FDG uptake in the posterior cingulate cortex could be relatively less affected

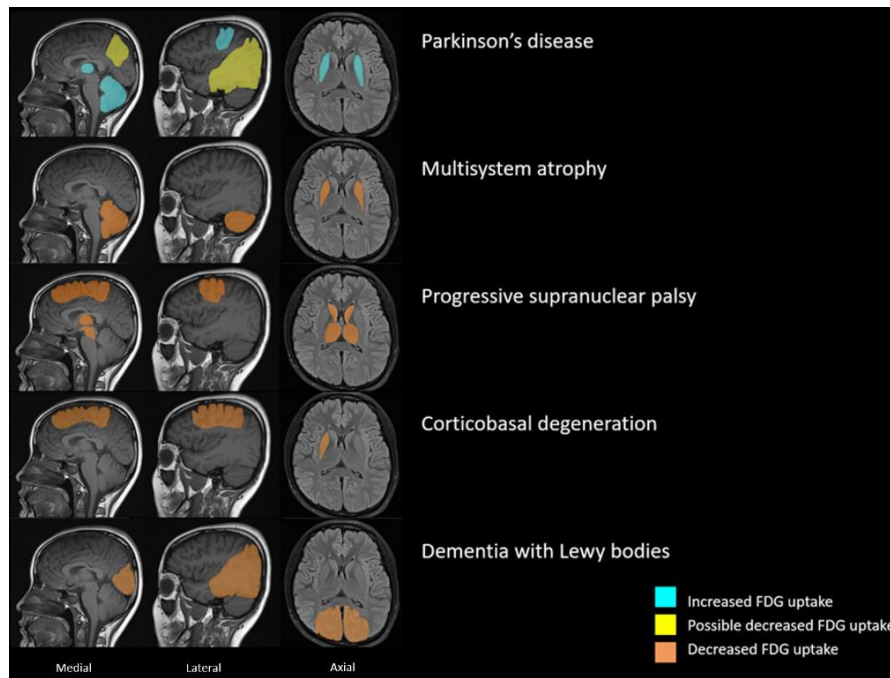


Figure 6. Summarization of brain glucose metabolic change in movement disorders

FDOPA PET

General Consideration

FDOPA or 6-fluoro-(¹⁸F)-L-3, 4-dihydroxy-phenyl -alanine is the amino acid analog precursor for dopamine. Therefore, FDOPA can detect abnormal dopaminergic presynaptic activity in movement disorder, confirming the clinical diagnosis of degenerative Parkinsonism (54).

Image Data Display

- The serial images would be re-oriented into the AC-PC plane and displayed in the axial plane.

The image orientation of symmetrical adjustment should be performed carefully to provide appropriate relative visual analysis to compare the degree of radiopharmaceutical uptake at both sides of the regional brain.

- Generally, the slice selection for image scaling can show the most prominent crescent shape of the putamen and the rounded shape of caudate.

Image Interpretation

The imaging findings and interpretation are shown in Table 19 and Figure 7.

Table 19. FDOPA PET Findings in Movement Disorder

Interpretation	Findings	Possible disease
Normal scan	the crescent shape of the putamen and rounded shape of caudate	- Normal
		- ET
	bilaterally symmetric and uniform thickness in both caudate and putamen (comma- or crescent-shaped)	- drug-induced parkinsonism
		- vascular parkinsonism
		- psychogenic parkinsonism
		- SWEDDs
Abnormal scan	reduction in size or thickness and shape of putamen or both putamen and caudate (unilaterally or bilaterally)	- PD
		- PDD
	Early stage: decrease at posterior putamen, contralateral to the affected side of the body	- PSP
	Late stage: Decrease at caudate and putamen	- CBD
		- MSA

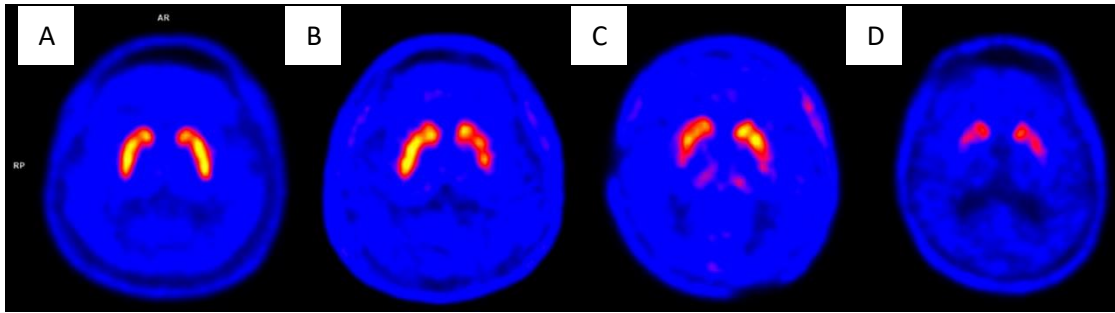


Figure 7. (A) A Normal FDOPA PET scan with intense symmetrical activity in bilateral striata. There is various severity of abnormal FDOPA PET scan in movement disorder as seen (B) showing asymmetrical decreased FDOPA uptake in left putamen and caudate with normal on the right side, (C) showing asymmetrical decreased FDOPA uptake in bilateral putamen and right caudate, and (D) showing decreased FDOPA uptake in bilateral caudate and putamen (more on the right side)

Differential Diagnosis in Parkinsonism

FDOPA PET imaging is not recommended to distinguish between PD and atypical Parkinsonism and among atypical Parkinsonism.

However, it can be used to differentiate between degenerative Parkinsonism and non-parkinsonian disorders, including essential tremor (ET) diagnosis, vascular Parkinsonism, drug-induced

Parkinsonism, and other non-parkinsonian disorders (55). Alteration of FDG uptake pattern in the striatum and the cortical brain area can be used to differentiate among degenerative Parkinsonism.

Summarization of nuclear medicine investigation and typical findings in movement disorder are shown in Figure 8.

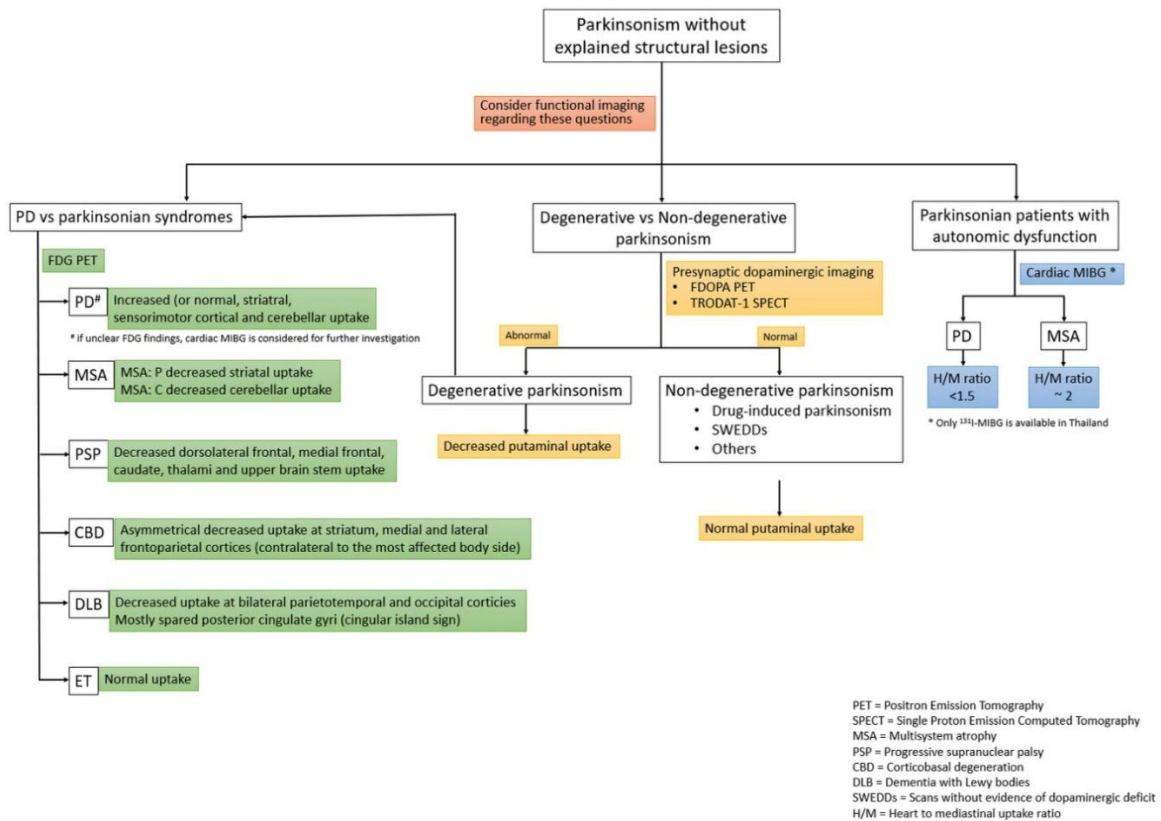


Figure 8. Summarization of nuclear medicine investigation and typical findings in movement Disorder

Reporting Format

Specific identification of each patient, referring clinician, date and time of the study, and the reporting physician’s signature must be provided in the report as part of quality assurance.

The recommended structural body of the report includes five major portions as follows:

1. History: The relevant history should be noted, such as specific symptoms and predominant side of movement disorder, the time duration

from the disease onset, current medication, and withdrawal time interval of drugs interfering with the study.

2. Indication

3. Technique

a. Radiopharmaceutical: type, dosage, and route of administration

b. Detailed imaging and processing techniques should be mentioned, including imaging quality and limitations. If certain specific software or anatomical co-registration is utilized, it should be additionally detailed.

c. Ancillary drugs (if used), e.g., type, route, and time of sedative drugs.

d. Serum glucose level for FDG PET.

4. Findings.

a. Quality of image, e.g., good, fair, poor (e.g., low count rate, motion).

b. Visual analysis, e.g., quality of uptake in caudate and anterior/posterior putamen of each side for TRODAT-1 or ¹⁸F-FDOPA. A term of hypometabolism and hypermetabolism

should be used for ¹⁸F-FDG. The location, extension, and severity should be reported.

c. Semi-quantitative analysis, e.g., type of normal database comparison used with Z-score cut-off, H/M ratio, and WOR for MIBG, TRODAT-1 SUR of caudate, putamen, and striatum, as well as the percentage of asymmetry should be reported.

d. Correlative imaging findings, if available, e.g., MRI, CT should be mentioned.

5. Interpretations/Impressions/Conclusions: The final impression from imaging data should be interpreted along with available clinical and correlative imaging data. If they seem discordant, direct discussion with the referring physician is recommended. Also, suggest additional study in some instances to clarify the diagnosis of the current study.

Pitfall and Error

Many factors, including biological and technical factors, affect image interpretation, as shown in Table 20(26, 28, 29, 44, 56-67).

Table 20. Pitfalls and Errors

Presynaptic dopaminergic imaging, including TRODAT-1 brain SPECT and FDOPA brain PET	
Biological factors	Technical factors
1. Normal variations <ul style="list-style-type: none"> - Decreased striatal binding in normal aging - Decreased striatal binding in postmenopausal women - Mild asymmetrical striatal uptake 2. Altered location or shape of the striatal structures due to anatomic lesions 3. Medications or drugs that may alter tracer binding (see Table XX in patient preparation section).	Technical factors <ol style="list-style-type: none"> 1. Artifacts <ul style="list-style-type: none"> - Patient movement - Misalignment causing artificial asymmetry 2. Variability of semi-quantitative analysis <ul style="list-style-type: none"> - Interobserver variability with manual ROI - Correction effects, e.g., attenuation, scatter, and partial volume effect 3. Suboptimal image display plane
Cardiac MIBG Imaging	
Biological factors	Technical factors
1. The progressive decline of MIBG uptake with aging 2. Age dependent metabolic and cardiovascular conditions 3. Heart diseases, e.g., ischemic heart disease, heart failure, dilated cardiomyopathy, myocardial damage, cardiac rhythm disease 4. Autonomic neuropathy or ganglionopathy, e.g., DM, post-polio syndrome 5. Medications <ul style="list-style-type: none"> - False positive (Amiodarone) - False-negative (Competitive compounds for NE transporter binding, calcium channel blocker, long-term amiodarone) 	1. Artifacts <ul style="list-style-type: none"> - Patient movement 2. Abnormal anatomy causing the erroneous position of mediastinal ROI
FDG PET	
Biological factors	Technical factors
1. Unexpected brain activity from external stimuli 2. Drug interference, e.g., psychotropic drugs, corticosteroid 3. The effects from sedation at the time of injection 4. Anatomical variations 5. Generalized reduced brain FDG uptake secondary to high blood glucose level 6. Recent radio- or chemotherapy 7. Age effect altered brain metabolism	1. Artifacts <ul style="list-style-type: none"> - Patient movement - Inappropriate processing, e.g., the reconstruction method 2. Non-continuous color table display 3. Level of contrast and background subtraction 4. Suboptimal image display plane, e.g., not true AC-PC plane, asymmetrical display 5. Partial volume effect on the corrected image

Statement and declaration

This guideline has got financial support from Global Medical Solution (GMS) Thailand, Siemens Healthcare, and Premier Business Inter.

All authors declare that they do not have any commercial or financial relationship that could be construed as the potential conflict of interest.

Conflict of interest

None.

Author's contributions

Conceptualization, S.T., B.Kh., T.K.; literature review, all authors; rate recommendation for indication in nuclear medicine imaging, P.L., O.P.; nuclear medicine imaging procedure, radiopharmaceuticals and interpretation, S.T., B.Kh., T.K., W.C., N.W., Ch.C., P.K., P.P., T.S., N.P., S.A.; pitfalls and errors, T.Th.; original draft preparation, T.K., B.Kh., S.T.; gathering external peer review, T.K.; writing review and editing, T.K., S.T.. All authors have read and agreed to the published version of this manuscript.

References

1. Tysnes OB, Storstein A. Epidemiology of Parkinson's disease. *J Neural Transm (Vienna)*. 2017; 124(8):901-5.
2. Eckert T, Barnes A, Dhawan V, Frucht S, Gordon MF, Feigin AS, et al. FDG PET in the differential diagnosis of parkinsonian disorders. *NeuroImage*. 2005; 26(3):912-21.
3. Hughes AJ, Daniel SE, Kilford L, Lees AJ. Accuracy of clinical diagnosis of idiopathic Parkinson's disease: a clinico-pathological study of 100 cases. *J Neurol Neurosurg Psychiatry*. 1992; 55(3):181-4.
4. Rajput AH, Rozdilsky B, Rajput A. Accuracy of clinical diagnosis in Parkinsonism--a prospective study. *Can J Neurol Sci*. 1991; 18(3):275-8.
5. Grimes D, Fitzpatrick M, Gordon J, Miyasaki J, Fon EA, Schlossmacher M, et al. Canadian guideline for Parkinson disease. *Cmaj*. 2019; 191(36):E989-e1004.
6. De la Fuente-Fernández R. Role of DaTSCAN and clinical diagnosis in Parkinson disease. *Neurology*. 2012; 78(10):696-701.
7. Postuma RB, Berg D, Stern M, Poewe W, Olanow CW, Oertel W, et al. MDS clinical diagnostic criteria for Parkinson's disease. *Mov Disord*. 2015; 30(12):1591-601.
8. Heinzl S, Berg D, Gasser T, Chen H, Yao C, Postuma RB, et al. Update of the MDS research criteria for prodromal Parkinson's disease. *Mov Disord*. 2019; 34(10):1464-70.
9. Godau J, Hussl A, Lolekha P, Stoessel AJ, Seppi K. Neuroimaging: current role in detecting pre-motor Parkinson's disease. *Mov Disord*. 2012; 27(5):634-43.
10. Stoessel AJ, Martin WW, McKeown MJ, Sossi V. Advances in imaging in Parkinson's disease. *Lancet Neurol*. 2011; 10(11):987-1001.
11. Politis M. Neuroimaging in Parkinson disease: from research setting to clinical practice. *Nat Rev Neurol*. 2014; 10(12):708-22.
12. Berardelli A, Wenning GK, Antonini A, Berg D, Bloem BR, Bonifati V, et al. EFNS/MDS-ES/ENS [corrected] recommendations for the diagnosis of Parkinson's disease. *Eur J Neurol*. 2013; 20(1):16-34.
13. Antonini A, Leenders KL, Vontobel P, Maguire RP, Missimer J, Psylla M, et al. Complementary PET studies of striatal neuronal function in the differential diagnosis between multiple system atrophy and Parkinson's disease. *Brain*. 1997; 120 (Pt 12):2187-95.
14. Brooks DJ, Salmon EP, Mathias CJ, Quinn N, Leenders KL, Bannister R, et al. The relationship between locomotor disability, autonomic dysfunction, and the integrity of the striatal dopaminergic system in patients with multiple system atrophy, pure autonomic failure, and Parkinson's disease, studied with PET. *Brain*. 1990; 113 (Pt 5): 1539-52.
15. Burn DJ, Sawle GV, Brooks DJ. Differential diagnosis of Parkinson's disease, multiple system atrophy, and Steele-Richardson-Olszewski syndrome: discriminant analysis of striatal ¹⁸F-dopa PET data. *J Neurol Neurosurg Psychiatry*. 1994; 57(3):278-84.
16. Arbizu J, Luquin MR, Abella J, de la Fuente-Fernandez R, Fernandez-Torron R, Garcia-Solis D, et al. [Functional neuroimaging in the diagnosis of patients with Parkinsonism: Update and recommendations for clinical use]. *Rev Esp Med Nucl Imagen Mol*. 2014; 33(4):215-26.
17. Postuma RB, Berg D, Stern M, Poewe W, Olanow CW, Oertel W, et al. MDS clinical diagnostic criteria for Parkinson's disease. *Movement Disorders*. 2015; 30(12):1591-601.
18. Goetz CG, Emre M, Dubois B. Parkinson's disease dementia: definitions, guidelines, and research perspectives in diagnosis. *Ann Neurol*. 2008; 64 Suppl 2:S81-92.
19. Brooks DJ. Diagnosis and management of atypical parkinsonian syndromes. *J Neurol Neurosurg Psychiatry*. 2002; 72 Suppl 1:I10-I6.
20. Respondek G, Stamelou M, Kurz C, Ferguson LW, Rajput A, Chiu WZ, et al. The phenotypic spectrum of progressive supranuclear palsy: a retrospective multicenter study of 100 definite cases. *Mov Disord*. 2014; 29(14): 1758-66.
21. Gilman S, Wenning GK, Low PA, Brooks DJ,

- Mathias CJ, Trojanowski JQ, et al. Second consensus statement on the diagnosis of multiple system atrophy. *Neurology*. 2008; 71(9):670-6.
22. McKeith IG, Boeve BF, Dickson DW, Halliday G, Taylor JP, Weintraub D, et al. Diagnosis and management of dementia with Lewy bodies: Fourth consensus report of the DLB Consortium. *Neurology*. 2017; 89(1):88-100.
 23. Montastruc JL, Llau ME, Rascol O, Senard JM. Drug-induced Parkinsonism: a review. *Fundam Clin Pharmacol*. 1994; 8(4):293-306.
 24. Kagi G, Bhatia KP, Tolosa E. The role of DAT-SPECT in movement disorders. *J Neurol Neurosurg Psychiatry*. 2010; 81(1):5-12.
 25. Schwingschuh P, Ruge D, Edwards MJ, Terranova C, Katschnig P, Carrillo F, et al. Distinguishing SWEDDs patients with asymmetric resting tremor from Parkinson's disease: a clinical and electrophysiological study. *Mov Disord*. 2010; 25(5):560-9.
 26. Djang DS, Janssen MJ, Bohnen N, Booij J, Henderson TA, Herholz K, et al. SNM practice guideline for dopamine transporter imaging with ¹²³I-ioflupane SPECT 1.0. *Journal of Nuclear Medicine*. 2012; 53(1):154-63.
 27. Bombardieri E, Giammarile F, Aktolun C, Baum RP, Bischof Delaloye A, Maffioli L, et al. ¹³¹I/¹²³I-metaiodobenzylguanidine (MIBG) scintigraphy: procedure guidelines for tumour imaging. *European journal of nuclear medicine and molecular imaging*. 2010; 37(12):2436-46.
 28. Varrone A, Asenbaum S, Vander Borght T, Booij J, Nobili F, Nagren K, et al. EANM procedure guidelines for PET brain imaging using [¹⁸F]FDG, version 2. *European journal of nuclear medicine and molecular imaging*. 2009; 36(12):2103-10.
 29. Waxman AD HK, Herholz K, Lewis DH, Herscovitch P, Minoshima S, Mountz JM, et al. Society of Nuclear Medicine Procedure Guideline for FDG PET Brain Imaging Version 1.0. 2009.
 30. Guedj E, Varrone A, Boellaard R, Albert NL, Barthel H, van Berckel B, et al. EANM procedure guidelines for brain PET imaging using [(18)F]FDG, version 3. *European journal of nuclear medicine and molecular imaging*. 2022; 49(2):632-51.
 31. Koopmans K, Ziekenhuis M, Eshuis GS. ¹⁸F DOPA PET/CT in Neuroendocrine Tumours and in Presynaptic Dopaminergic Deficits.
 32. Chondrogianis S, Marzola MC, Rubello D. (18)F-DOPA PET/computed tomography imaging. *PET clinics*. 2014; 9(3):307-21.
 33. Yang T, Wang L, Li Y, Cheng M, Jiao J, Wang Q, et al. (131)I-MIBG myocardial scintigraphy for differentiation of Parkinson's disease from multiple system atrophy or essential tremor in Chinese population. *J Neurol Sci*. 2017; 373:48-51.
 34. Jokinen P, Helenius H, Rauhala E, Bruck A, Eskola O, Rinne JO. Simple ratio analysis of ¹⁸F-fluorodopa uptake in striatal subregions separates patients with early Parkinson disease from healthy controls. *Journal of nuclear medicine: official publication, Society of Nuclear Medicine*. 2009; 50(6): 893-9.
 35. Sun Y, Chen Z, Liu C, Tang J, Wang Y, Zhang G, et al. Phase I clinical study with different doses of (99m)Tc-TRODAT-1 in healthy adults. *Ann Nucl Med*. 2020; 34(3):212-9.
 36. Xu D, Zhu W, Huo L, Zhu S, Li F, Wang H. Validation of Iodine-131-meta-iodobenzyl-guanidine cardiac scintigraphy in Parkinsonism: A preliminary study. *Parkinsonism Relat Disord*. 2018; 50:69-73.
 37. Ziessman MD HA OMM, Janis P, Thrall MD, James H. *Nuclear Medicine The Requisites 4th Edition* Saunders 2013. 2013.
 38. Yeh SH, Tsai CJ, Yu TH, Chiang YH, Lin SZ, Peng NJ, et al. ^{99m}Tc-TRODAT-1 SPECT Revealed That Striatal Dopamine Transport Availability Significantly Decreases in Late Mid-Aged Healthy Taiwanese and Then Remains Stable. *Clin Nucl Med*. 2022; 47(3): 201-8.
 39. Treglia G, Cason E, Stefanelli A, Cocciolillo F, Di Giuda D, Fagioli G, et al. MIBG scintigraphy in differential diagnosis of Parkinsonism: a meta-analysis. *Clinical autonomic research: official journal of the Clinical Autonomic Research Society*. 2012; 22(1):43-55.
 40. Orimo S, Suzuki M, Inaba A, Mizusawa H. ¹²³I-MIBG myocardial scintigraphy for differentiating Parkinson's disease from other neurodegenerative parkinsonism: a systematic review and meta-analysis. *Parkinsonism Relat Disord*. 2012; 18(5): 494-500.
 41. Yoshita M. Differentiation of idiopathic Parkinson's disease from striatonigral degeneration and progressive supranuclear palsy using iodine-123 meta-iodobenzyl-guanidine myocardial scintigraphy. *J Neurol Sci*. 1998; 155(1):60-7.
 42. Druschky A, Hilz MJ, Platsch G, Radespiel-Troger M, Druschky K, Kuwert T, et al. Differentiation of Parkinson's disease and multiple system atrophy in early disease stages by means of I-123-MIBG-SPECT. *J Neurol Sci*. 2000; 175(1):3-12.
 43. Jindahra P, Vejajiva A, Witoonpanich R, Sirisriro R, Sritara C, Pulkes T. Differentiation of dementia with Lewy bodies, Alzheimer's disease and vascular

- dementia by cardiac ^{131}I -meta-iodobenzyl-guanidine (MIBG) uptake (preliminary report). *Journal of the Medical Association of Thailand = Chotmaihet thangphaet*. 2004; 87(10):1176-81.
44. Lavalaye J, Booij J, Reneman L, Habraken JB, van Royen EA. Effect of age and gender on dopamine transporter imaging with [^{123}I] FP-CIT SPET in healthy volunteers. *European journal of nuclear medicine*. 2000; 27(7):867-9.
 45. Yang BH GY, Wu TH. Tc-99m TRODAT-1 SPECT imaging protocol guideline. *Annals of Nuclear Medicine and Molecular Imaging* 2016; 29:45-53.
 46. Yin TK, Lee BF, Yang YK, Chiu NT. Differences of various region-of-interest methods for measuring dopamine transporter availability using $^{99\text{m}}\text{Tc}$ -TRODAT-1 SPECT. *Scientific World Journal*. 2014; 2014: 837439.
 47. Huang WS, Lin SZ, Lin JC, Wey SP, Ting G, Liu RS. Evaluation of early-stage Parkinson's disease with $^{99\text{m}}\text{Tc}$ -TRODAT-1 imaging. *Journal of nuclear medicine: official publication, Society of Nuclear Medicine*. 2001; 42(9): 1303-8.
 48. Berti V, Pupi A, Mosconi L. PET/CT in diagnosis of movement disorders. *Annals of the New York Academy of Sciences*. 2011; 1228: 93-108.
 49. Brajkovic L, Kostic V, Sobic-Saranovic D, Stefanova E, Jecmenica-Lukic M, Jesic A, et al. The utility of FDG-PET in the differential diagnosis of Parkinsonism. *Neurological research*. 2017; 39(8):675-84.
 50. Hellwig S, Amtage F, Kreft A, Buchert R, Winz OH, Vach W, et al. [(18)F]FDG-PET is superior to [(123)I]IBZM-SPECT for the differential diagnosis of parkinsonism. *Neurology*. 2012; 79(13):1314-22.
 51. Meyer PT, Frings L, Rucker G, Hellwig S. (18)F-FDG PET in Parkinsonism: Differential Diagnosis and Evaluation of Cognitive Impairment. *Journal of nuclear medicine: official publication, Society of Nuclear Medicine*. 2017; 58(12):1888-98.
 52. Song IU, Ha SW, Yang YS, Chung YA. Differences in Regional Glucose Metabolism of the Brain Measured with F-18-FDG-PET in Patients with Essential Tremor According to Their Response to Beta-Blockers. *Korean journal of radiology*. 2015; 16(5):967-72.
 53. Walker Z, Gandolfo F, Orini S, Garibotto V, Agosta F, Arbizu J, et al. Clinical utility of FDG PET in Parkinson's disease and atypical Parkinsonism associated with dementia. *European journal of nuclear medicine and molecular imaging*. 2018; 45(9):1534-45.
 54. Chevalme Y-M, Montravers F, Vuillez J-P, Zanca M, Fallais C, Oustrin J, et al. FDOPA- (^{18}F): a PET radiopharmaceutical recently registered for diagnostic use in countries of the European Union. *Brazilian Archives of Biology and Technology*. 2007; 50:77-90.
 55. Highlight of prescribing information F-DOPA [Internet]. 2019. Available from: https://www.accessdata.fda.gov/drugsatfda_docs/label/2019/200655s000lbl.pdf.
 56. Booij J, Kemp P. Dopamine transporter imaging with [^{123}I] FP-CIT SPECT: potential effects of drugs. *European journal of nuclear medicine and molecular imaging*. 2008; 35(2):424-38.
 57. Tatsch K, Någren K, Nobili F, Walker Z, Darcourt J, Booij J, Tatsch K, et al. EANM procedure guidelines for brain neuro-transmission SPECT using (123) I-labelled dopamine transporter ligands, version 2. 2009.
 58. Badiavas K, Molyvda E, Iakovou I, Tsolaki M, Psarrakos K, Karatzas N. SPECT imaging evaluation in movement disorders: far beyond visual assessment. *European journal of nuclear medicine and molecular imaging*. 2011; 38(4):764-73.
 59. Linke R, Gostomzyk J, Hahn K, Tatsch K. [^{123}I] IPT binding to the presynaptic dopamine transporter: variation of intra- and inter observer data evaluation in parkinsonian patients and controls. *European journal of nuclear medicine*. 2000; 27(12):1809-12.
 60. Morton RJ, Guy MJ, Clauss R, Hinton PJ, Marshall CA, Clarke EA. Comparison of different methods of DatSCAN quantification. *Nuclear medicine communications*. 2005; 26(12):1139-46.
 61. Tissingh G, Booij J, Bergmans P, Winogrodzka A, Janssen AG, van Royen EA, et al. Iodine-123-N-omega-fluoropropyl-2beta-carbomethoxy-3beta-(4-iodophenyl) tropane SPECT in healthy controls and early-stage, drug-naive Parkinson's disease. *Journal of nuclear medicine: Official publication, Society of Nuclear Medicine*. 1998; 39(7):1143-8.
 62. Jacobson AF, Travin MI. Impact of medications on MIBG uptake, with specific attention to the heart: comprehensive review of the literature. *Journal of Nuclear Cardiology*. 2015; 22(5):980-93.
 63. Nakajima K, Yoshita M, Matsuo S, Taki J, Kinuya S. I-123 MIBG sympathetic imaging in Lewy-body diseases and related movement disorders. *QJ Nucl Med*. 2008; 52: 378-87.
 64. Skowronek C, Zange L, Lipp A. Cardiac ^{123}I -MIBG Scintigraphy in Neurodegenerative Parkinson Syndromes: Performance and Pitfalls in Clinical Practice. *Front Neurol*. 2019; 10: 152.
 65. Tsujikawa K, Hasegawa Y, Yokoi S, Yasui K,

- Nanbu I, Yanagi T, et al. Chronological changes of ^{123}I -MIBG myocardial scintigraphy and clinical features of Parkinson's disease. *Journal of Neurology, Neurosurgery & Psychiatry*. 2015; 86(9):945-51.
66. Borsatti E. PET and PET/CT-a clinical guide, EC Lin, A. Alavi (Eds.), Thieme, New York/Stuttgart (2009), 312 pages, 505 illustrations, US \$79.95, ISBN: 978-1-60406-153-6. Elsevier; 2010.
67. Tepmongkol S, Srikijvilaikul T, Vasavid P. Factors affecting bilateral temporal lobe hypometabolism on ^{18}F -FDG PET brain scan in unilateral medial temporal lobe epilepsy. *Epilepsy & Behavior*. 2013; 29(2): 386-9.